

Tetramethylfluoroformamidinium hexafluorophosphate (TFFH)/ benzyltriphenylphosphonium dihydrogen trifluoride (PTF): a unique reagent for the conversion of carboxylic acids to the corresponding alcohols as well as hydroxamic acids

Ayman El-Faham,^{a*} Sherine N. Khattab,^a and Mohamed Abdul-Ghani^b

^aDepartment of Chemistry, Faculty of Science, Alexandria University,
Ibrahimia 21321, P.O. Box 426, Alexandria, EGYPT

^bDepartment of Chemistry, Faculty of Science, Beirut Arab University,
P.O. Box 11-5020, Beirut, LEBANON
E-mail: Aymanel_faham@hotmail.com

Abstract

The present work describes the utility of TFFH (**1**) for the conversion of carboxylic acids to the corresponding alcohols **4** as well as the hydroxamic acids **5** in the presence or absence of the fluoride additive PTF (**3**). The addition of the fluoride additive **3** during the activation step improves the yield and purity of the products, due to the fast activation to the acid fluoride derivatives as observed from IR examination.

Keywords: Tetramethylfluoroformamidinium hexafluorophosphate (TFFH), benzyltriphenylphosphonium dihydrogen trifluoride (PTF), alcohols, hydroxamic acids

Introduction

Halouronium salts are frequently described as good reagents especially for coupling of sterically hindered *N*-substituted or *N,N*-disubstituted amino acids, which is not easily carried out under standard conditions.¹ Tetramethylfluoroformamidinium hexafluorophosphate **1** (TFFH), was introduced by Carpino and El-Faham² as a benign substitute for the corrosive cyanuric fluoride **2**. Cyanuric fluoride **2** was reported as a mild reagent, suitable for the preparation of acyl fluorides even when unsaturated double bonds, hydroxyl groups, or aromatic rings are present.^{3,4} TFFH has been shown to act as a coupling reagent that proceeds *via in situ* conversion to an acid fluoride.^{2,5} Acyl fluorides, having a greater stability than the corresponding chlorides towards neutral oxygen nucleophiles such as water and methanol, appear to be of equal reactivity toward anionic nucleophiles and amines.⁶ Recently, Fmoc-,⁷ Boc-, or *Z*-amino acid fluorides⁸ have been found to be stable and fast acylating reagents for peptide bond formation. Moreover, TFFH

Table 1. Reduction of carboxylic acids and *N*-protected amino acids to alcohols

RCOOH	RCH ₂ OH ^a	Yield (%) ^b TFFH (TFFH/PT F)	mp. °C (<i>lit.</i> , °C)	[α] _D ²⁵ ^c
4a C ₆ H ₅ -COOH	5a C ₆ H ₅ CH ₂ OH	78 (93)	—	—
4b 2-NO ₂ -C ₆ H ₄ -COOH	5b 2-NO ₂ -C ₆ H ₄ CH ₂ OH	80 (94)	—	—
4c C ₆ H ₅ -CH ₂ -COOH	5c C ₆ H ₅ CH ₂ -CH ₂ OH	79 (91)	—	—
4d 9-Fluorene carboxylic acid	5d 9-Fluorenylmethanol	75 (92)	97 (100) ^{3a}	—
4e Boc-Phe-COOH	5e Boc-Phe-CH ₂ OH	84(89)	94-95 (94.5) ^{3b}	-23.5 (c 1, CHCl ₃)
4f Z-Phe-COOH	5f Z-Phe-CH ₂ OH	80 (84)	91-92 (90-91) ^{3b}	-41.1 (c 1, MeOH)
4g Fmoc-Leu-COOH	5g Fmoc-Leu-CH ₂ OH	88 (91)	130-131 (130) ^{3b}	-17.4 (c 1, CHCl ₃)
4h Fmoc-Glu(OBn)-COOH	5h Fmoc-Glu(OBn)CH ₂ OH	81(85)	91 (91-92) ^{3b}	-9.0 (c 1, CHCl ₃)

^a The products **5a-h** were identified by their physical constants, IR, ¹H NMR data and comparison with data reported in the literature.


^b Yields are of isolated and purified products

^c The present method (TFFH alone or TFFH/PTF) proceeds with retention of optical purity as indicated by comparison of the specific rotation values with those reported in the literature.³

N-Protected amino alcohols are important synthetic intermediates, particularly in the synthesis of peptides that are potent inhibitors of proteases,¹² and are also incorporated at the C-terminal of biological active peptides.¹³ The conversion of amino acids into alcohols by chemoselective reduction of their corresponding mixed anhydrides¹⁴ and *N*-carboxyanhydrides¹⁵ with sodium borohydride has also been reported. In the present work, a variety of *N*-protected amino acids **4e-h** were converted rapidly into alcohols **5e-h** by *in situ* reduction of their corresponding acid fluorides with sodium borohydride in high yield. There was not much difference between TFFH alone or the TFFH/PTF method (Table 1), because the symmetric anhydride was not observed during IR examination.

The present work also describes a useful method for the preparation of a variety of hydroxamic acids using TFFH (**1**) alone, or in the presence of the fluoride additive PTF (**3**). Carboxylic acids **6a-b** as well as *N*-protected amino acids **6c-i** were activated for 5 min with TFFH alone, or the TFFH/PTF mixture, in the presence of triethylamine (Et₃N) or diisopropylethyl amine (DIEA) as a base in dichloromethane (DCM) followed by addition of

hydroxylamine.HCl at 0 °C. The reaction mixture was warmed to room temperature and then stirred at room temperature overnight. After workup, the products were obtained in high yields and purity (Table 2, Scheme 2). There was no difference in the yield and purity for urethane *N*-protected amino acids with the use of TFFH alone or TFFH/PTF. However, there was a difference in the yield in case of aryl carboxylic acids; TFFH/PTF gave better yields and purity than TFFH alone (Table 2).

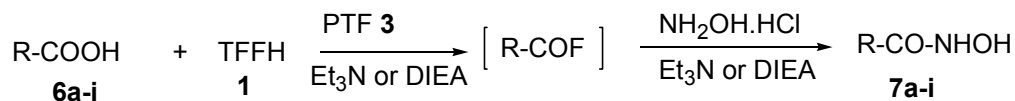

Scheme 2

Table 2. Preparation of aryl hydroxamic acids and *N*-protected amino hydroxamic acids using TFFH or TFFH/PTF*

R-COOH	R-NHOH ^a	Yield (%) ^b TFFH (TFFH/PTF)	mp. (°C) (<i>lit.</i> , °C)	[α] _D ²⁵ ^c
6a C ₆ H ₅ -COOH	7a C ₆ H ₅ CONHOH	81 (95)	123 (123-124) ¹⁷	
6b C ₆ H ₅ -CH ₂ -COOH	7b C ₆ H ₅ -CH ₂ -NHOH	80 (91)	Oil (oil) ¹⁷	
6c Boc-Leu-OH	7c Boc-Leu-NHOH	90(94)	Oil (oil) ¹⁷	-8.9 (c 1, MeOH)
6d Boc-Phe-OH	7d Boc-Phe-NHOH	89(93)	Oil (oil) ¹⁷	+25.7 (c 1, EtOH)
6e Boc-Glu(OBn)-OH	7e Boc-Glu(OBn)-NHOH	88(91)	83 (82-83) ¹⁷	+27.8 (c 1, MeOH)
6f Z-Phe-OH	7f Z-Phe-NHOH	80 (84)	91-93 (90-91) ¹⁷	-38.6 (c 0.5, MeOH)
6g Z-Val-OH	7g Z-Val-NHOH	90 (92)	Oil (oil) ¹⁷	+10.5 (c 1, MeOH)
6h Fmoc-Phe-OH	7h Fmoc-Phe-NHOH	88 (91)	150-151 (150-152) ¹⁷	-39.3 (c 1, DMF)
6i Fmoc-Leu-OH	7i Fmoc-Leu-NHOH	94(96)	115 (114-116) ¹⁷	-25.4 (c 1, DMF)

*All compounds **7a-i** gave a positive test for hydroxamic acids on treatment with iron (III) chloride solution.¹⁶

^a IR and ¹H NMR spectra of compounds **7a-i** were in accordance with their structure

^b Yields are of isolated and purified products

^c The present method proceeds with retention of optical purity as indicated by comparison of the specific rotation values with those reported in the literature.¹⁷

Conclusions

The TFFH/PTF procedure provides a general, rapid, and convenient method for the conversion of carboxylic acids in the presence of NaBH₄ into alcohols as well as to hydroxamic acids. The present method proceeds also with retention of optical purity as indicated by comparison of the specific rotation values with those reported in the literature.

Experimental Section

General Procedures. Melting points was obtained in open capillary tubes using a melting heating apparatus and are uncorrected. NMR spectra were recorded using a Bruker 300 MHz or a Jeol 500 MHz instrument with TMS as internal standard. Infrared spectra were obtained using a Shimadzu 8300 series Fourier Transform instrument. Optical rotations were measured at 25 °C. All asymmetric amino acids derivatives were of the L-configuration and were purchased from Fluka Chemical Co. All solvents and chemical were of reagent grade and used without further purification. Silica gel 60 (70-230 mesh, Merck) was used for column chromatography.

General Procedure for the Preparation of Alcohols (5a-h): TFFH or TFFH/PTF (2 mmol, 0.528 g/0.824 g respectively) was added to a stirred solution of the acid **4a-d** (2 mmol) and Et₃N (4 mmol, 0.404 g) or *N*-protected amino acid **4e-h** (2 mmol) and diisopropylamine (DIEA, 4 mmol, 0.516 g) in CH₂Cl₂ (5 mL), and kept under a N₂ atmosphere at 0 °C. After the mixture was stirred at 0 °C for about 1 h, ice-cold water (10 mL) was added along with CH₂Cl₂ (20 mL). The organic layer was separated, and the aqueous layer was extracted once with CH₂Cl₂ (10 mL). The combined organic layer was washed with ice-cold water (10 mL), dried (MgSO₄), and concentrated under reduced pressure to a small volume (2 mL). NaBH₄ (76 mg, 2 mmol) was added in one portion, and MeOH (2 mL) was then added dropwise over a period of 15-20 min at rt. The reaction mixture was neutralized with 1N H₂SO₄, and the organic solvents were evaporated under pressure. The residue was treated with EtOAc (10 mL) and H₂O (5 mL); the organic layer was separated, and the aqueous layer was extracted with EtOAc (2 x 10 mL). The combined organic layers were washed consecutively with 1N H₂SO₄ (5 mL) and H₂O (2 x 10 mL) and dried (MgSO₄), the solvent was evaporated under reduced pressure. The residue was purified by distillation or column chromatography using EtOAc/hexane (1:1) as an eluent.

General Procedure for the Preparation of Hydroxamic Acids (7a-i): TFFH or TFFH/PTF (2 mmol, 0.528 g/0.824 g respectively) was added to a stirred solution of the acid **6a-b** (2 mmol) and Et₃N (4 mmol, 0.404 g) or *N*-protected amino acid **6c-i** (2 mmol) and diisopropylamine (DIEA, 4 mmol, 0.516 g) in CH₂Cl₂ (10mL) and was activated for 5-10 min at 0 °C. Hydroxylamine hydrochloride (140 mg, 2 mmol) and Et₃N or diisopropylamine (2 mmol) were added to the reaction mixture at 0 °C. After the addition, the reaction mixture was warmed to

room temperature, and then stirred overnight, filtered through Celite, and the organic phase washed with 1N HCl (3x5 mL), a saturated NaCl solution (2x5 mL), and dried (MgSO₄). The solvent was evaporated under reduced pressure to yield compounds **7a-i** that were isolated without further purification.

Acknowledgements

Professor L. A. Carpino, Professor of Organic Chemistry, University of Massachusetts, Amherst, U. S. A. is thanked for his support and advice.

References and Notes

[†]**Abbreviations used:** TFFH = Tetramethylfluoroformamidinium hexafluorophosphate; DCM = dichloromethane, EtOAc = ethylacetate; PTF = benzyl triphenylphosphonium dihydrogen trifluoride; DIEA = diisopropylethylamine; Et₃N = triethylamine; Z-Phe- = *N*-benzyloxycarbonyl phenylalanine; Z-Val- = *N*-benzyloxycarbonyl valine; Boc-Phe- = *N*-*tert*-butyloxycarbonyl phenylalanine; Boc-Leu- = *N*-*tert*-butyloxycarbonyl leucine; Fmoc-Leu- = 9-fluorenylmethyloxycarbonyl leucine. Amino acids are abbreviated and designated following the rules of the IUPAC-IUB Commission of Biochemical Nomenclature (*J. Biol. Chem.*, **1972**, 247, 977).

1. (a) Tung, R. H; Rich, D. H. *J. Am. Chem. Soc.* **1985**, *107*, 4342. (b) Coste, J.; Dufour, M.-N.; Pantaloni, A.; Castro, B. *Tetrahedron Lett.* **1990**, 669. (c) Frérot, E; Coste, J.; Pantaloni, A.; Dufour, M.-N. ; Jouin, P. *Tetrahedron* **1991**, *47*, 259. (d) Coste, J.; Frérot, E.; Jouin, P. *Tetrahedron Lett.* **1991**, *32*, 1967. (e) Wenschuh, H.; Beyermann, M.; Krause, E.; Carpino, L. A.; Bienert, M. *Tetrahedron Lett.* **1993**, 3733.
2. (a) Carpino, L. A.; El-Faham, A. *J. Am. Chem. Soc.* **1995**, *117*, 5401. (b) El-Faham, A. *Org. Prep. Proced. Int.* **1998**, *30*, 477. (c) El-Faham, A. *Chemistry Lett.* **1998**, 671.
3. Olah, G.A.; Nojima, M.; Kerekes, I. *Synthesis* **1973**, 487.
4. (a) Carpino, L. A. *J. Org. Chem.* **1980**, *45*, 4250. (b) Bertho, J.; Loffet, A.; Pinel, C.; Reuther, F.; Sennyey, G. *Tetrahedron Lett.* **1991**, *32*, 1303. (c) Kokotos, G.; Noula, C. *J. Org. Chem.* **1996**, *61*, 6994.
5. (a) Wenschuh, H. ; Beyermann, M.; Krause, E.; Brudel, M.; Winter, R.; Schümann, M.; Carpino, L. A.; Bienert, M. *J. Org. Chem.* **1994**, *59*, 3275. (b) Carpino, L. A.; Beyermann, M.; Wenschuh, H.; Bienert, M. *Acc. Chem. Res.* **1996**, *29*, 268.
6. (a) Bender, M. L.; Jones, M. L. *J. Org. Chem.* **1962**, *27*, 3771. (b) Swai, C. G.; Scott, C. B. *J. Am. Chem. Soc.* **1953**, *75*, 246.

7. Carpino, L. A.; Sadat-Aalae, D.; Chao, H. G.; DeSelms, R. H. *J. Am. Chem. Soc.* **1990**, *112*, 9651.
8. (a) Carpino, L. A.; Mansour, E. S. E.; Sadat-Aalae, D. *J. Org. Chem.* **1991**, *56*, 2614. (b) Carpino, L. A.; Masnour, E. S. E.; El-Faham, A. *J. Org. Chem.* **1993**, *58*, 4162.
9. Boss, U.; Pedersen, B.; Christensen, J. B. *Synth. Commun.* **1998**, *28*, 1223.
10. El-Faham, A.; Abdul-Ghani, M. *Org. Prep. Proced. Int.* **2003**, *35*, 369.
11. Carpino, L. A.; Ionescu, D.; El-Faham, A.; Beyermann, M.; Henklein, P.; Hanay, C.; Wenschuh, H.; Bienert, M. *Org. Lett.* **2003**, *5*, 975.
12. (a) Thompson, R. *Methods Enzymol.* **1977**, *46*, 220. (b) Sharma, R. P.; Gore, M. G.; Akhtar, M. *Chem. Commun.* **1979**, 875. (c) Bakker, A. V.; Jung, S.; Spencer, R. W.; Vinick, F.; Faraci, W. S. *Biochem. J.* **1990**, *271*, 559.
13. Pless, J.; Bauer, W.; Cardinaux, F.; Closse, A.; Hauser, D.; Huguenia, R.; Roemer, D.; Buescher, H. H.; Hill, R. C. *Helv. Chim. Acta.* **1979**, *62*, 398.
14. (a) Kokotos, G. *Synthesis*, **1990**, 299. (b) Rodriguez, M.; Linares, M.; Doulut, S.; Heitz, A.; Martinez, J. *Tetrahedron Lett.* **1991**, *32*, 923.
15. Fehrentz, J. A.; Califano, J. C.; Amblard, M.; Loffet, A.; Martinez, J. *Tetrahedron Lett.* **1994**, *35*, 569.
16. *Vogel's Textbook of Practical Organic Chemistry* Longman: London 5th Edn., 1989, p 1222.
17. Giacomelli, G.; Porcheddu, A.; Salari, M. *Org. Lett.* **2003**, *5*, 2715.